


Rear wheel size:
 Tire type/size:
 Caster wheel size:
 Caster fork:

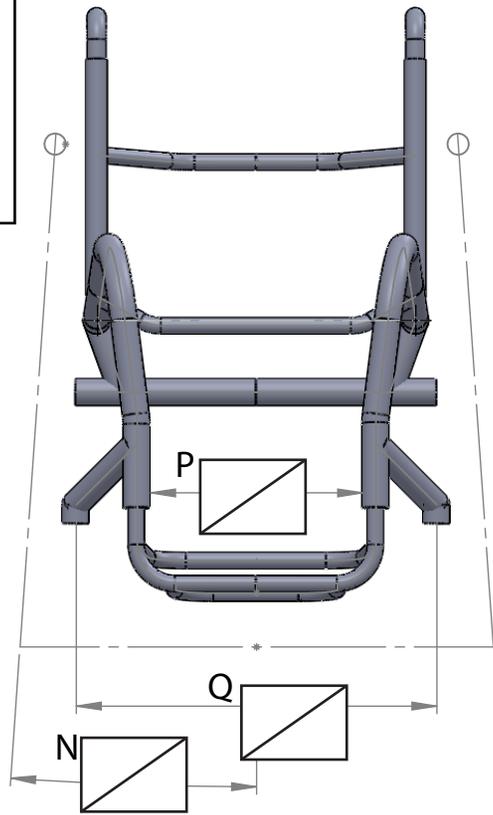

Measurements V & W
 are new, not in
 measurement video

OLD
 NEW

	Measurement Form
Name:	_____
Date:	_____
Phone:	_____
Email:	_____

- A- Measure from ground to top of back upholstery or rigid backrest
- B- With both rear wheels against wall, measure from wall to back of back post at the bottom of the post where it meets the seat tube/frame rails
- C- Measure from the wall to the center of your rear wheel. This is easiest to measure to the center of the quick release axle.
- D- Measure the overall length of chair by placing the end of the tape measure against the wall and measuring to the front of the footrest
- E- Measure from the wall to the center of your caster housings
- F- Measure the depth of your seat pan. Measure from back of seat pan to front of seat pan
- G- Measure height of back posts from the seat. Measure from the bottom of the back posts to the top
- H- Measure the rear seat height. Measure from the ground to the back of the seat sling
- I- Measure the front seat height. Measure from the ground to the front of the seat sling
- J- Measure the height of the footrest. Measure from the ground to the top of the footplate at the front of the footrest
- K- Measure the angle of the footrest using an angle measuring app on your cell phone
- L- Measure the angle of the front of your wheelchair frame using an angle measuring app on your cell phone
- M- Measure the angle of the back post using an angle measuring app on your cell phone. If you have a Rigid back, please also measure the angle of the backrest itself
- N- Measure the camber angle of your wheels using an angle measuring app on your cell phone
- O- Measure seat width. Measure from the inside of one clothing guard to the inside of the other. If you do not have clothing guards, measure outside of frame tube to outside of frame tube
- P- Measure the width of the front of your frame from inside of tube to inside of tube where the footrest inserts into the main frame
- Q- Measure the width of your casters from center of caster housing to center of caster housing
- R- Measure width of footrest
- S- Measure depth of footrest from the back of your footrest to the front
- T- If you have ergo seating, measure the flat in your seat rails from back of your seat sling to the point the seat rails start to bend
- U- If you have push canes, measure from the ground to the top of your push canes
- V- If you have seat taper, measure from the outside of left seat tube to outside of right seat tube at the front of your seat sling.
- W- If you are looking to get a double foot hoop on your new chair, please measure the distance you want the heel bar back from the front of the footrest